

Medical beams



Why choose MZ Liberec?



Design - Ergonomics - Maintenance

- Light aluminium construction, smooth surface simplifies product cleaning and disinfection
- Colour may be adapted to the interior design
- Broad range of additional accessories available




Functions and features

MZ Liberec medical beams are unmistakable for their design, variability and modularity. A custom-made design allows you to adapt the appearance to the overall interior and to the needs of the medical staff.

- A. electrical channel
- B. gas channel
- C. rails for slider
- D. lighting

Technical properties

- Medical beams are primarily designed for specialized healthcare facilities (ICU, OR).
- They contain outlets for medical gases, high and low-voltage electrical sockets.
- Beams can be equipped with a wide range of accessories thanks to the rail system.
- Multi-bed options.

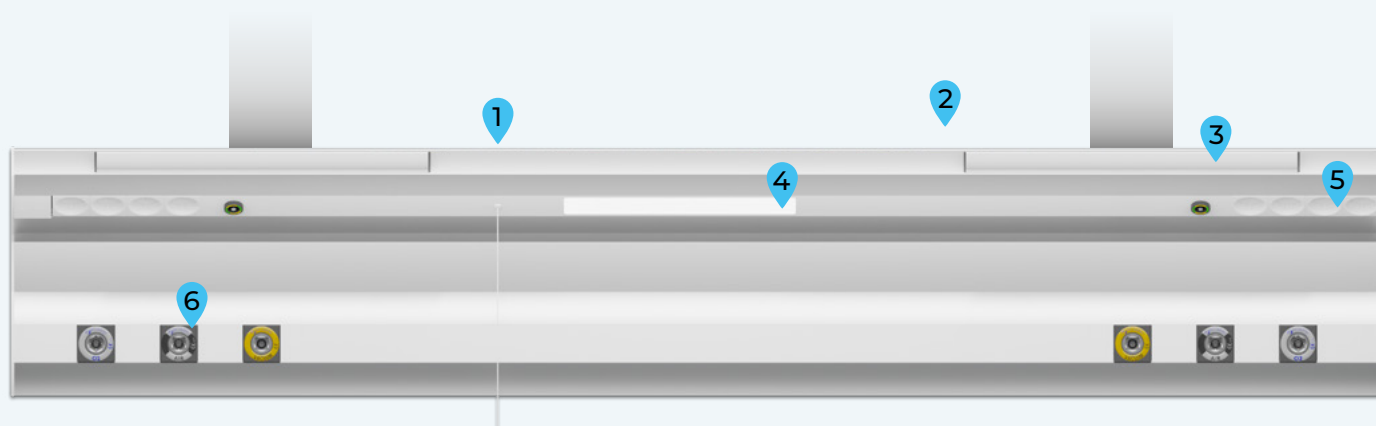


energy-saving LED lighting

- 1. indirect lighting
- 2. night lighting
- 3. medical rail

- 4. direct lighting
- 5. electric socket
- 6. gas outlets


Integrated electrical equipment

The embedded design makes cleaning and disinfection simple

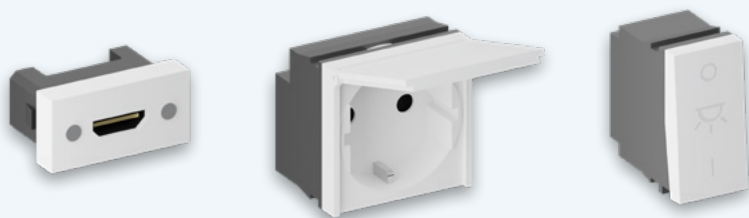


Gas outlets

Standards: ČSN, DIN, AFNOR, BS, AGA + other

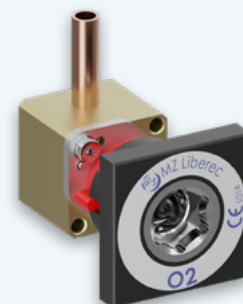

Models - ZMP07

ZMP07 - single sided

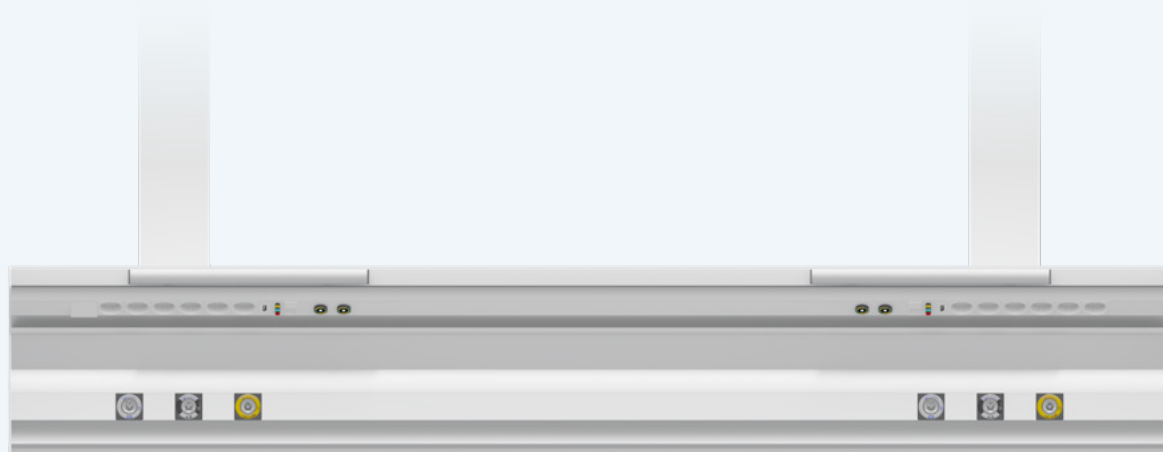

ZMP07 - double sided




ZMP07 - laminar



Technical parameters - ZMP07

Models ZMP	ZMP07 - single sided	ZMP07 - double sided	ZMP - laminar
Standard length (pro 1 bed) *multi-bed design available	1950 mm	1950 mm	3000 mm (one side)
Length of supporting columns	2000 mm		
Weight (standard length)	max. 150 kg	max. 175 kg	max. 300 kg
Lighting (LED)	direct / indirect / night		
Electrical devices	Outlets in different colours with or without light, designed according to national standards, RJ45 data connector, phone and nurse call connection, etc.		
Medical gas outlets	All types of terminal units according to national standards (DIN, AFNOR, BS, AGA...) outlets can led from both sides at the back		
Material design	Steel construction of the supporting profile with powder coating/ Anodized aluminum profiles with removable covers that can be in RAL colour		

Possible colour design

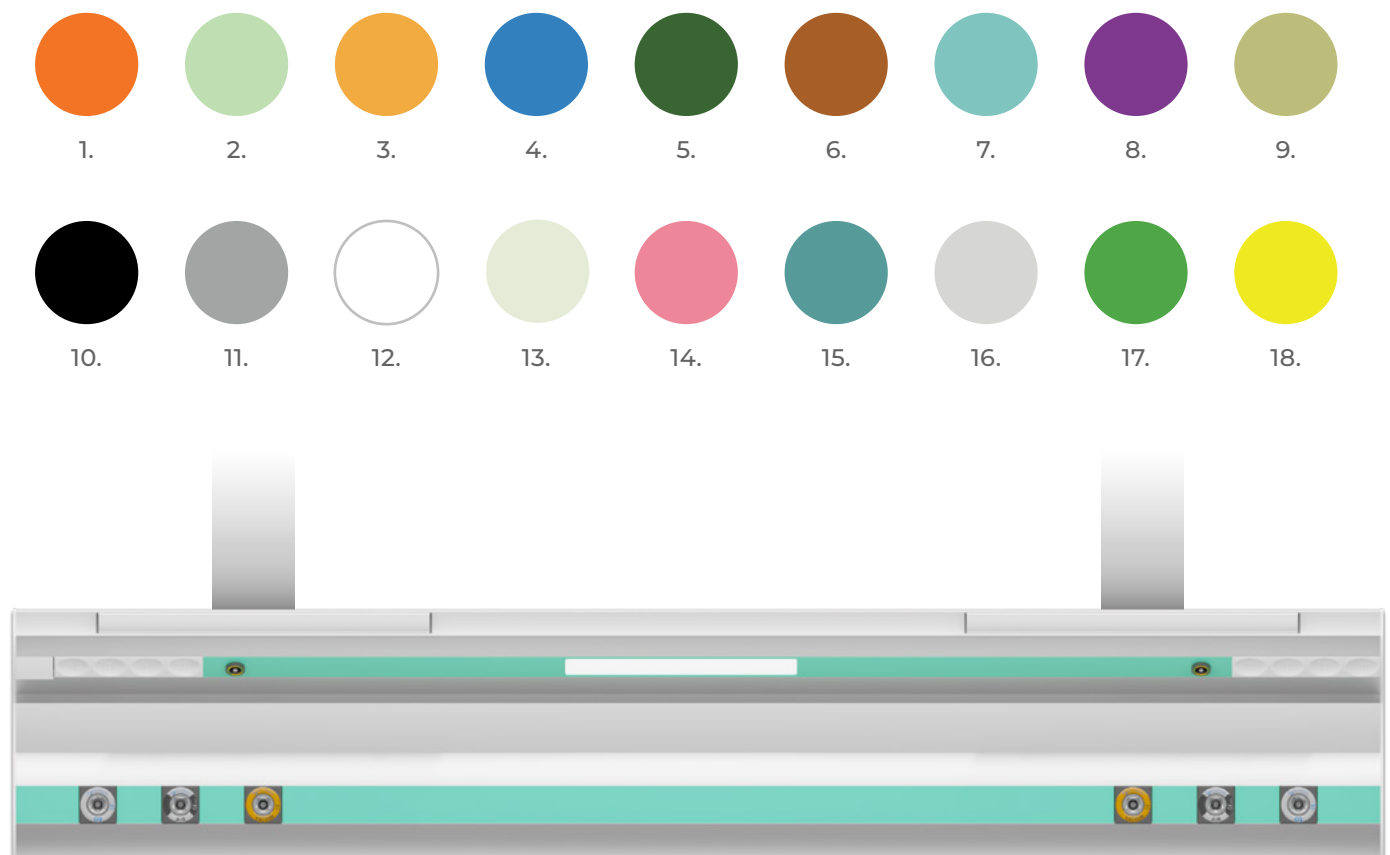

 - We adapt the colour design to your interior

 - Antibacterial surface

Accessories


Shelves

Bar mounted shelves serve as a storage space for common medical supplies and aids so they remain close at hand. At the same time, the arm movement control button may be positioned on a shelf.




Shelf wire baskets

Provides a storage space for surgical instruments and other medical supplies during surgeries and for patients in ICU, ARD units.

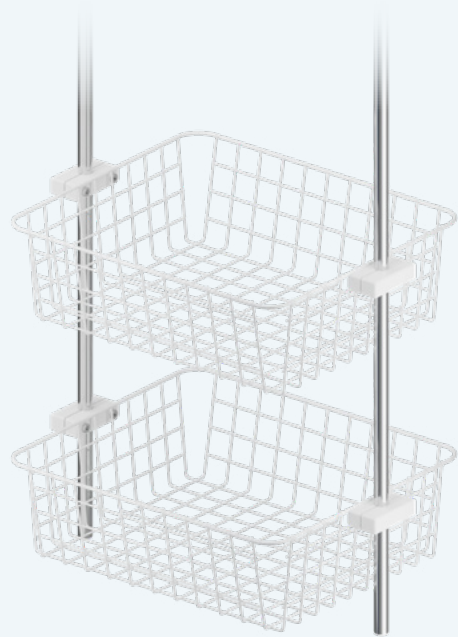

Support arms for additional shelves and infusions


Many accessories can be placed on the hangers, especially infusion racks, stands, infusion pump holders and monitors.

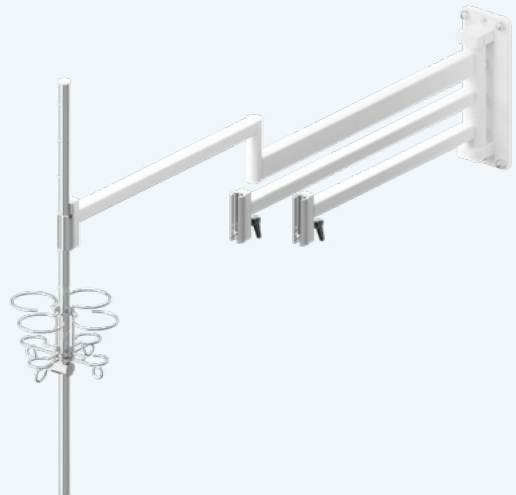

Accessories

SPR double compact arms with electric terminal

These are a huge help wherever the patient's health condition requires the use of many devices at once. The advantage is the additional electrical power they provide.

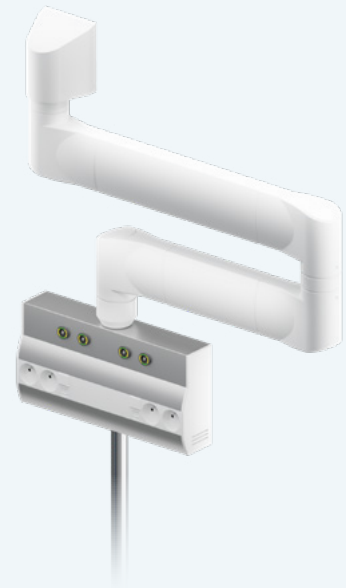


Rail slider

Bar mounted shelves serve as a storage space for common medical supplies and aids so they remain close at hand. At the same time, the arm movement control button may be positioned on a shelf.



Support arm (single / double)

Monitors, infusion racks and other accessories can be installed on the hanger.

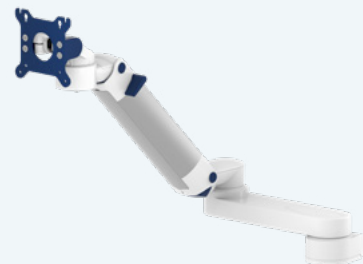

LED light

The LED light is designed for reading or simple examinations conducted at the patient's bed. Its ergonomic design is the idea of daily aid for patients and medical staff.



Additional accessories are available in our catalogue and on our website.



www.mzliberec.cz

Rudník 4, 543 72 Rudník u Vrchlabí, CZ

ID:47306581, VAT: CZ47306581

Copyright © MZ Liberec, a.s., 2023

Verze: 2311

