

MEDICAL CATALOGUE

www.rimsa.it

MEDICAL CATALOGUE

INDEX

Our History	4
Main Features	6
Unica	8
U29	10
N Series	12
E Series	14
Quattroluci LED	16
Pentaled 30E Light	17
Pentaled 12/28	18
Saturno LED	20
Prima LED	21
Alfa LED	22
L88 LED TW	23

OUR HISTORY

TRADITION FOR INNOVATION

RIMSA, established in 1936 by Palmino Longoni, was initially a mechanical workshop dedicated to repairing typewriters and the like; hence, the acronym R.I.M.S.A. (Riparazione di Macchine da Scrivere e Affini).

The transition from repair workshop to production facility took place in the 1940s, when Mr. Palmino decided to give shape to a product of his own. Since then, RIMSA dedicated itself to the design and development of balanced arm lamps. Company growth resulted in an expansion of the product range with the introduction of magnifying and fluorescent lamps. Starting in the post-war period, RIMSA began making a name for itself in the electronics, goldsmithery, dentistry and industrial sectors.

In the 80s, RIMSA began focusing closely on the surgical lighting sector and, in April 1983, the Milan Trade Fair Authority awarded RIMSA the prize for the design of a halogen surgical lamp. Research in the medical field continued and in March 1992 the Milan Chamber of Commerce awarded the company the prestigious "Technological Innovation" qualification certificate for the design of the Stellare open-spoke surgical lamp for operating theatres with laminar-flow.

In 2002, RIMSA developed the world's first LED operating light, at a time when this technology was still in its infancy.

In 2017, RIMSA patented "UNICA" the world's first surgical light with no glare.

RIMSA Società Benefit is a certified ISO 9001, 13485, 14001 and 27001, B Corp. RIMSA belongs to the WHP (workplace health promotion) project. All our products are marked CE.

1936
Establishment

1943
Palmino Longoni becomes a member of the Italian Inventors Guild

1945
RIMSA's first balanced arm lamp

1956
Historical Company certificate

1971
Purposedly designed CNC and lathe light

1983
RIMSA's first surgical light

1992
Technological Innovation Award

2002
RIMSA pioneers the first LED surgical light

2017
Unica is the first LED surgical light with no glare

2024
RIMSA certified B Corp

DETERMINATION AND PASSION: **THE ART OF INNOVATING**

RIMSA IS A HISTORICAL COMPANY,
NOT AN OLD ONE.

Backed by its history, traditions and pride, RIMSA has always put “Uniqueness” at the centre of its organization, focusing on the promotion of human resources, technological updating, straightforward management, and product quality. Together, all these elements lead to the achievement of the corporate purpose, identified as follows: continuity and development of the Company, professional growth and staff development, research and innovation, and acquisition of new markets.




MAIN FEATURES


LIGHT EMITTING SURFACE

The main characteristic of scialytic, or surgical, lamps can be found in their name. The term scialytic originally comes from the Greek for 'devoid of shadow': scialytic lamps are therefore shadowless. Over the decades many different lighting technologies have been adopted. In 2002, at the "Medica" trade fair in Düsseldorf, RIMSA presented the first LED powered scialytic lamp in the world. The LED technology guarantees a lower energy consumption and a longer service life compared to any other lighting source. On the other hand, since LEDs are spaced between each other on the reflector, the light emitting surface generated by DIRECT LED light is smaller compared to that of halogen light.


In order to compensate for this problem, RIMSA studied INDIRECT LED light, which, together with the other advantages offered by LEDs, allowed LED technology to become the main light source in the medical field.

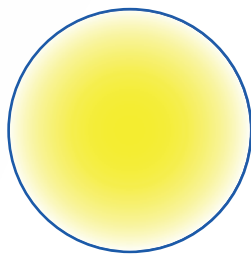


Halogen indirect light
100%

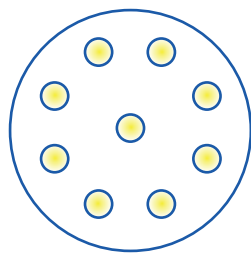

LED direct light
60%

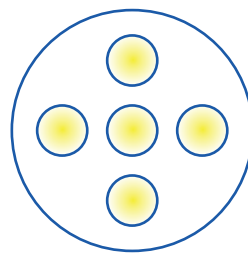

LED indirect light
80%

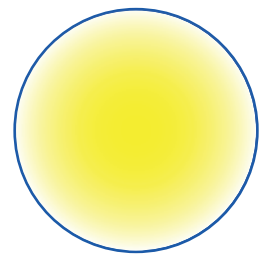

2R double reflection (UNICA)
100%


% of light emitting surface of the cupola depending on the lighting technology used

INDIRECT LIGHT

Conscious of the disadvantages caused by direct illumination, RIMSA designed and implemented a solution using INDIRECT light for its own products: the beam of light that is produced by the diode is intercepted by a parabola that reflects the beams, merging them into one. Therefore, indirect light allows for a greater illuminated surface area whilst using fewer LEDs.

Added benefits of indirect light also include a lower luminous degradation, a greater scialytic effect than that of a direct light solution, and a lower temperature generated on the printed circuit board, guaranteeing a longer service life.


Further developing the advantages of indirect light, in 2017 RIMSA patented its 2R technology (Double reflection technology) which guarantees the maximisation of the light emitting surface area.

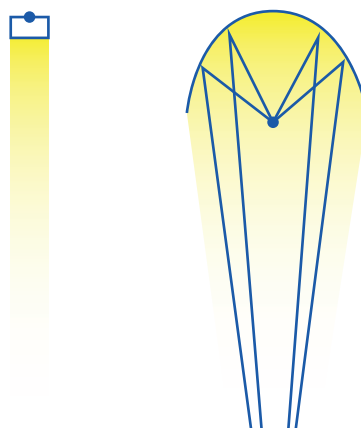

Visualisation of the principles of direct light (left) and indirect light (right)

2R: DOUBLE REFLECTION TECHNOLOGY

The presence of LEDs along two circumferences positioned at the extremity of two modules of mirrors creates a double reflection illumination capable of obtaining a light emitting surface area that is incomparable to any other type of reflection. With the patented 2R technology, the light emitting surface area matches the area of the reflector.

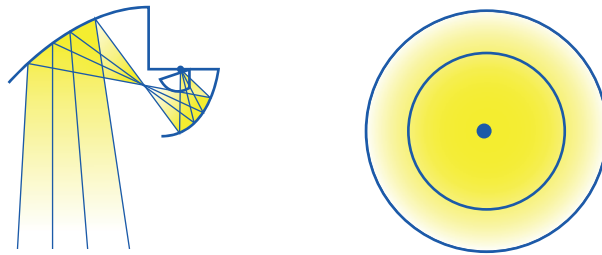

Representation of double reflection technology (left) and maximized light emitting surface (right)

GLARING EFFECT

Considering the duration required for some operations and the potential gravity of human error in these situations, it is essential to guarantee a light that doesn't stress, tire, or blind the surgeon, even temporarily. The most irritating effect that scalytic lamps can create is glare: the sensation of being blinded is particularly dangerous when we consider that these lamps are required to produce up to 160.000 lux.

The lamps of the Unica series, thanks to the 2R Technology, are capable of completely removing the sensation of glare caused by surgical illumination.

COMPENSATION AREA

An additional light source makes it possible to expand the lit field at the edges without affecting the light intensity at the centre (Ec). The compensation area is a ring of less intense light surrounding the operating field, which allows a progressive and gradual adaptation of the light in order to remove the glare between the central area of the surgical focus and the peripheral one.

DEEP 3D LIGHT

The centre of the lamp is fitted with an additional LED module specifically designed to cast deep light. This feature guarantees that the surgeon will be able to operate with perfect 3D lighting, especially in cavities.

UNICA

THE ONLY LAMP THAT COMPLETELY ELIMINATES GLARE

With the models of the Unica series, RIMSA has obtained the complete elimination of glare. Thanks to its double reflection technology, it is possible to obtain a complete suppression of shadows, the maximization of the light emitting surface and an absolute glare free light.

UNICA 860
UNICA 520



Verified EPD: UNICA provides concrete data on the entire life cycle, from raw material extraction to final disposal. EPD declaration no. 4791865733.101.1, verified by UL Solutions.





CEILING SINGLE
UNICA520SO



MOBILE
UNICA520PI



CEILING DOUBLE
UNICA520+520



CEILING SINGLE
UNICA860SO



CEILING DOUBLE
UNICA860+520

Performances

	860	520
Illumination (Ec) at 1 m distance	160.000 lux	160.000 lux
Light head diameter	86 cm	52 cm
Colour temperature (7 selections / continuous)	3.800 to 5.400 K	3.800 to 5.400 K
Light field diameter	170 - 400 mm	170 - 380 mm
Average LED life	> 60.000 hours	> 60.000 hours
Light intensity control	3 - 100 %	3 - 100 %

Main features

Double reflection: zero glare

Deep 3D light

Compensation area

Light always on focus

Future proof: all light heads supplied camera ready as a standard

Accessories




Double yoke




Camera:
4K cable
4K Wi-Fi
FULL-HD cable
FULL-HD Wi-Fi



Monitor

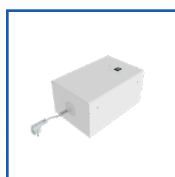

Battery group

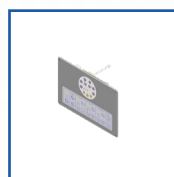

Wall control
panel



IR remote
control



Integration
system

U29

A COMFORTABLY AFFORDABLE SURGICAL LIGHT

U29 combines the indirect light technology consolidated by the success of the PentaLED Series with the elegant and compact design of Unica.

U29 is a high-performance lamp, strongly engineered and comfortably affordable, so that any surgeon can benefit from it.



Verified EPD: U29 provides concrete data on the entire life cycle, from raw material extraction to final disposal.

EPD declaration no. 4791865733.101.1, verified by UL Solutions.



CEILING SINGLE
U29SO



MOBILE
U29PI



CEILING DOUBLE
U29+29

Performances

U29

Illumination (Ec) at 1 m distance	160.000 lux
Light head diameter	52 cm
Colour temperature	4.500 – 5.000 K
Light field diameter	140 – 230 mm
Average LED life	> 60.000 hours
Light intensity control	5 – 100 %

Main features

- Indirect light
- Deep 3D light
- Compensation area
- Light always on focus
- Future proof: all light heads supplied as standard camera ready

Accessories




Double yoke




Camera:
4K cable
4K Wi-Fi
FULL-HD cable
FULL-HD Wi-Fi



Monitor

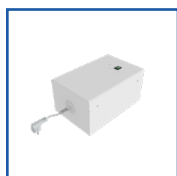

Battery group

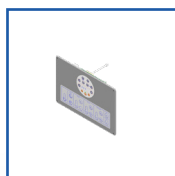

Wall control
panel



IR remote
control



Integration
system

PENTALED N SERIES

MANUAL DIAMETER ADJUSTMENT

The PentaLED N Series represented a breakthrough in the medical lighting field, being the first LED surgical light ever displayed to the public. The focus function is activated by the surgeon rotating the sterile central handle: this manual focalization inside the sterile area grants a precise and immediate control of the light field.

PENTALED 63N PENTALED 30N





CEILING SINGLE
PENTA30NSO



MOBILE
PENTA30NPI



WALL MOUNTED
PENTA30NPA



CEILING DOUBLE
PENTA30N+30N



CEILING SINGLE
PENTA63NSO



CEILING DOUBLE
PENTA63N+30N



CEILING DOUBLE
PENTA63N+63N

Performances

	63 N	30 N
Illumination (Ec) at 1 m distance	160.000 lux	160.000 lux
Light head diameter	63 cm	40 cm
Colour temperature	4.500 – 5.000 K	4.500 – 5.000 K
Diameter adjustment	Manual	Manual
Light field diameter	160 – 300 mm	140 – 280 mm
Average LED life	> 60.000 hours	> 60.000 hours
Light intensity control	5 – 100 %	5 – 100 %

Main features

- Indirect light
- Deep 3D light
- Manual diameter and focus adjustment

Accessories




Double yoke




PENTA63N ONLY
Camera:
FULL-HD cable
4K cable



Monitor

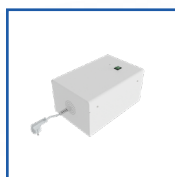

Battery group

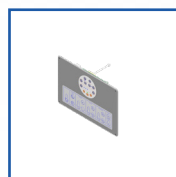

Wall control
panel



IR remote
control



Integration
system

PENTALED E SERIES

ELECTRONIC DIAMETER ADJUSTMENT

The electronic adjustment allows the operator to easily swap between 2 diameter sizes from the keyboard without modifying the light intensity at the center.

PENTALED 30E





CEILING SINGLE
PENTA30ESO



MOBILE
PENTA30EPI



WALL MOUNTED
PENTA30EPA



CEILING DOUBLE
PENTA30E+30E

Performances

30 E

Illumination (Ec) at 1 m distance	160.000 lux
Light head diameter	40 cm
Colour temperature	4.500 – 5.000 K
Light field diameter	140 – 260 mm
Average LED life	> 60.000 hours
Light intensity control	5 – 100 %

Main features

Indirect light
Deep 3D light
Compensation area
Light always on focus

Accessories




Double yoke




Camera:
FULL-HD cable
4K cable



Monitor

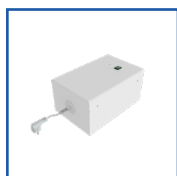

Battery group

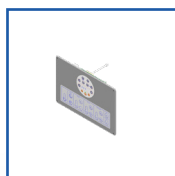

Wall control
panel



IR remote
control



Integration
system

QUATTROLUCI LED

IDEAL FOR SURGERY, UP TO 160.000 LUX
FROM A 60 CM REFLECTOR

This model is especially recommended for operating theatres where the surgeon needs a light and thin lamp to avoid interference with other overhead equipment. Excellent for oral and maxillofacial surgery, and plastic surgery. The newly updated multi-circular dome, compact and wieldy, consists of 36 LED lenses divided into 4 modules, which ensure shadow suppression and three-dimensional light.



CEILING SINGLE
QUATTROSO



MOBILE
QUATTROPI



WALL MOUNTED
QUATTROPA



CEILING DOUBLE
QUATTROSOX2

Performances


QUATTROLUCI

Illumination (Ec) at 1 m distance	160.000 lux
Light head diameter	60 cm
Colour temperature	4.500 K
Light field diameter	300 mm
Average LED life	> 60.000 hours
Light intensity control	15 – 100%

Main features

Deep 3D light
Direct light

Accessories

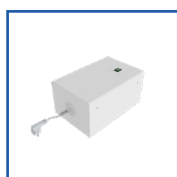

Battery group

PENTALED 30E LIGHT

QUALITY LIGHT WITH A COMPETITIVE PRICE

PentaLED 30E LIGHT combines the well-known cupola from PentaLED 30E with a lighter structure, making it even more cost-effective without compromising on light quality.



CEILING SINGLE
PENTA30ELSO



MOBILE
PENTA30ELPI



WALL MOUNTED
PENTA30ELPA



CEILING DOUBLE
PENTA30EL+30EL

Performances


30E LIGHT

Illumination (Ec) at 1 m distance	160.000 lux
Light head diameter	40 cm
Colour temperature	4.500 - 5.000 K
Light field diameter	140 - 260 mm
Average LED life	> 60.000 hours
Light intensity control	5 - 100 %

Main features

Indirect light
Deep 3D light
Compensation area
Light always on focus

Accessories

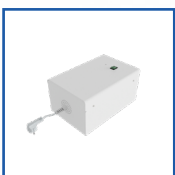

Battery group

PENTALED 12/28

UNPARALLELED PERFORMANCE TECHNOLOGY,
THE BEST LAMP FOR OUTPATIENT CLINICS AND MINOR SURGERY.

Their compact dimensions and extremely handy structures are permeated by the strong determination of RIMSA's dedication, a synonym of high technology, quality and performance since forever. The thin dome, with its two convenient side grips, ensures easier positioning and adjustment and reduces overall dimensions.

Choose the lamp that best meets your necessities: Pentaled 12 comes with indirect light; Pentaled 28 on the other hand adopts a direct light technology and has the possibility to manually adjust the light field diameter.



PENTALED 12

PENTALED 28



CEILING SINGLE
PENTA12SO
PENTA28SO



MOBILE
PENTA12PI
PENTA28PI



WALL MOUNTED
PENTA12PA
PENTA28PA



CEILING DOUBLE
PENTA12+12
PENTA28+28

Performances

12

28


Illumination (Ec) at 1 m distance	100.000 lux	120.000 lux
Light head diameter	40 cm	40 cm
Colour temperature	4.500 K	4.500 – 5.000 K
Diameter adjustment	Fixed	Manual
Light field diameter	160 mm	110 – 330 mm
Average LED life	> 60.000 hours	> 60.000 hours
Light intensity control	20 – 100 %	20 – 100 %

Main features

Indirect light – PENTALED 12

Direct light – PENTALED 28

Accessories

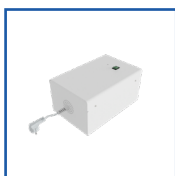

Battery group

SATURNO LED

QUALITY LIGHT WITH A SMALL REFLECTOR

This surgical-type lamp is suitable for minor non-invasive surgeries, for gynaecology and for the emergency room. The fact that the beams are close together (reflector size 195 mm), means that they do not have to be focused. The lamp is very easy to move thanks to the lightness of the aluminium support structure.



CEILING SINGLE
SATSON-LED



MOBILE
SATPIN-LED



WALL MOUNTED
SATPAN-LED



CEILING DOUBLE
SATSONX2-LED

Performances

	SATURNO
Illumination (Ec) at 1 m distance	60.000 lux
Light head diameter	19,5 cm
Colour temperature	4.000 - 4.500 K
Light field diameter	260 mm
Average LED life	> 60.000 hours
Light intensity control	20 - 100%

Main features

Deep 3D light
Direct light

Accessories




Battery group

OBSERVA SERIES

PRIMA LED

The first and only examination lamp with the possibility to choose an ambient light in addition to traditional concentrated light.

The Ergo-Spring balancing system makes PRIMALED very stable and easy to operate. PRIMALED is ideal for any type of installation, from the outpatient department to the intensive care unit.

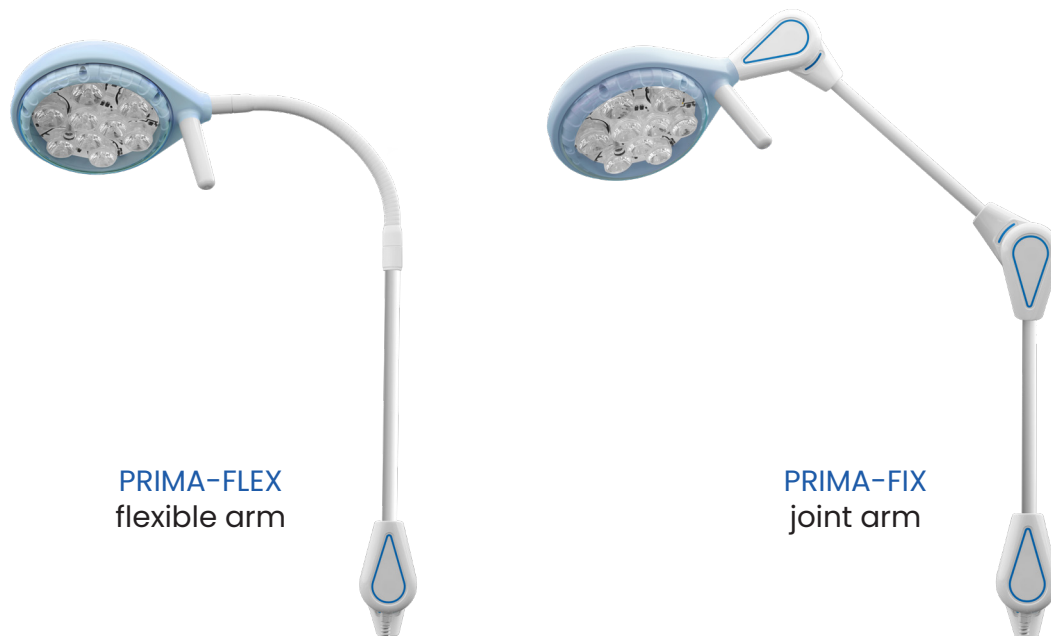

PRIMA-FLEX
flexible arm

PRIMA-FIX
joint arm

Performances

Illumination at 0,50 m distance	115.000 lux
Light head diameter	19,5 cm
Colour temperature	4.000 - 4.500K
Light field diameter	150 mm
Average LED life	> 60.000 hours
Light intensity control	20 - 100 %

PRIMA

Fixing systems



S11
Table clamp



S12MED
Wall clamp



RLBI
Mobile base



Z400819
Rail bar clamp

Accessories



Battery group

OBSERVA SERIES

ALFA LED

Three LED light sources with coinciding lenses and borosilicate glass protection shield ensure intense, deep light with low power consumption and minimal heat radiation. Each LED is equipped with a resistor to ensure continuous use of the lamp even in the rare event of an LED failing. Thanks to these three independent light sources, shadows can be reduced.



ALFA-FLEX
flexible arm



ALFA-FIX
joint arm

Performances

Illumination at 0,50 m distance	70.000 lux
Light head diameter	9,5 cm
Colour temperature	4.000 K
Light field diameter	140 mm
Average LED life	> 60.000 hours
Light intensity control	4 - 100 %

ALFA

Fixing systems



S11
Table clamp



S12MED
Wall clamp



RLBI
Mobile base for
ALFA-FIX



RLALFA
Mobile base for
ALFA-FLEX



Z400819
Rail bar clamp




Battery group

Accessories

OBSERVA SERIES

L88 LED TW

Especially suitable for dermatological use and wherever magnifying in general is needed, this model features a biconvex magnifying lens in optical glass with Ø 120 mm. It features a polycarbonate screen for protecting the light source. The lamp is equipped with a 3 dioptries lens; 5 dioptries lens available on request. The lamp is now available in the tunable white version: when the lamp is switched off, colour temperature can be adjusted from 2.700 K to 6.000 K. With the lamp switched on, light intensity can be adjusted.

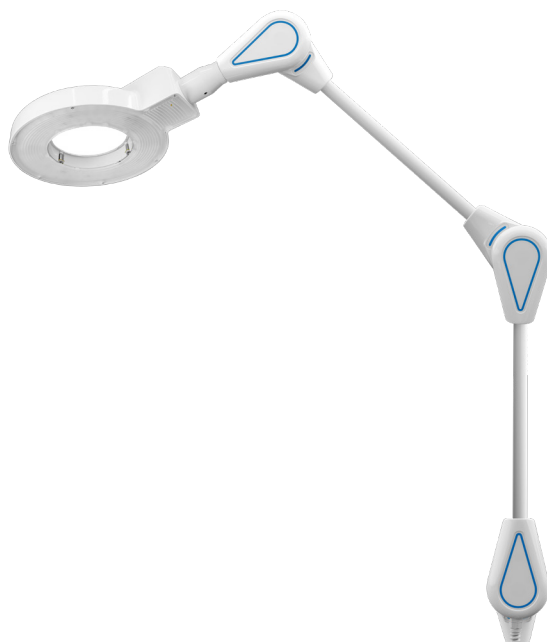

Performances

L88

Illumination at 0,50 m distance	2.250 lux
Light head diameter	23 cm
Colour temperature	2.700 to 6.000 K
Light field diameter	200 mm
Average LED life	> 60.000 hours
Light intensity control	5 - 100 %

Fixing systems



S11
Table clamp



S12MED
Wall clamp



RLBI
Mobile base



Z400819
Rail bar clamp

Accessories



Battery group



Rimsa - Brightening ideas

Via Monterosa 18 / 22 - 20831 Seregno (MB) - Italy
Tel. (+ 39) **0362 325709** | E-mail **info@rimsa.it**

RIMSA retains the right to improve the products in the catalogue without notice. Reproduction in part or in whole is forbidden.

www.rimsa.it



Research & Components
HANDMADE IN ITALY

